AT52336

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En

Thank you for purchasing our product. Please read this instruction manual carefully before use. Do not use the product before reading and fully understanding this manual. If you do not understand the warnings, notes, and recommendations, contact a healthcare professional or the seller to avoid product damage or personal injury. **WARNING:** All parts must be checked for damage. Before first use, thoroughly inspect the technical condition of the product and ensure that all parts are complete and properly assembled (if applicable, pay particular attention to screws, nuts, rivets and other fastening and structural components). Make sure that the product is not defective or damaged. If any damage or missing parts are found, DO NOT use the product. For further information, contact the seller.

INTENDED PATIENT GROUP

Individuals suffering from diseases, dysfunctions, or injuries for the treatment, rehabilitation, or compensation of which this product is intended (see the section on the intended use of the product in this user manual).

The product may be purchased independently by the user or based on the recommendation of a physician, therapist, or other specialist. Whether the product is purchased independently or recommended by a physician/therapist/other specialist, available sizes, required functions and variants of the product, indications and contraindications for use, as well as the information provided by the manufacturer must be taken into account.

INTENDED USE: An electric wheelchair is a vehicle structurally designed for the mobility of a person with a disability. It is powered by a motor whose design limits the maximum speed to walking pace.

INDICATIONS

Severe mobility limitations: Individuals with significant limitations in muscle strength, for example due to neurological or muscular diseases or pain syndromes, which prevent or significantly hinder independent mobility.

Neurodegenerative diseases: Individuals suffering from conditions such as multiple sclerosis, Parkinson's disease, muscular dystrophies, or ALS, which lead to progressive loss of muscle control and strength.

Fatigue: Individuals for whom prolonged walking or manual operation of a wheelchair is excessively demanding due to illness or health condition, causing rapid fatigue and exhaustion.

Age-related limitations: Older people who, due to age-related reductions in mobility and strength, may benefit from an electric wheelchair to maintain independence and mobility.

Rehabilitation and recovery periods: Individuals undergoing rehabilitation after injuries or surgeries who are temporarily unable to move independently may use an electric wheelchair as part of the treatment and recovery process.

CONTRAINDICATIONS

Physical or mental limitations (e.g. visual impairment) that prevent safe handling of the product.

Individuals who have difficulty judging distances or spatial relationships may have problems safely maneuvering an electric wheelchair.


Electric wheelchairs require a certain level of manual ability and coordination to operate the joystick or other control devices. Persons who are unable to operate these controls effectively may not be able to use the electric wheelchair safely.

Although many individuals with cognitive impairments can safely use electric wheelchairs, in some cases severe impairments may prevent safe operation. Individuals who are unable to maintain a stable seated position may require additional support systems to safely use an electric wheelchair. Epilepsy or other conditions that may lead to sudden loss of consciousness or reaction may constitute a contraindication. If a person has impaired motor coordination or is unable to independently operate the wheelchair controls, this may be a contraindication. Individuals with disorders of consciousness that may affect their ability to safely operate an electric wheelchair may also be unsuitable for its use. In some cases, serious heart conditions or blood pressure disorders may limit the safe use of an electric wheelchair, especially if emotions or stress may trigger adverse reactions. Conditions such as advanced multiple sclerosis or other diseases affecting muscle control may require specialized adaptive solutions in wheelchairs. Problems with depth perception, spatial perception, or general orientation may significantly reduce the safety of using an electric wheelchair. Generally poor health or recent surgeries, especially involving the spine or lower limbs, may require postponing the use of the wheelchair until health has improved.

LIST OF COMPONENTS

Wheelchair, battery, charger, remote control, user manual

- 1.Backrest
- 2.Controller
- 3.Footrests
- 4.Front wheel
- 5.Push handles
- 6.Armrest
- 7.Seat
- 8.Frame
- 9.Battery
- 10.Drive wheel motor
- 11.Charger

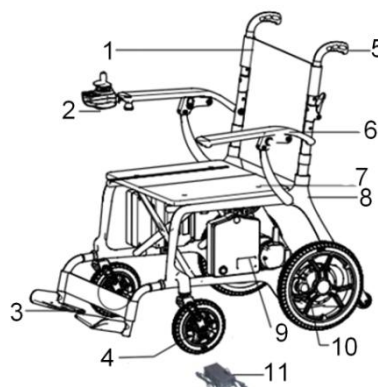TECHNICAL SPECIFICATIONS

Material: Full carbon fiber frame

Dimensions (unfolded): Length 103 cm × Width 72.5 cm × Height 102 cm

Dimensions (folded): Length 72.5 cm × Width 48 cm × Height 78.5 cm

Motor: 2 × 24 V / 250 W brushless motors

Controller: 24 V / 35 A LED controller

Charger: 24 V / 2 A

Battery: 24 V / 13 Ah lithium-ion battery

Brake: Electromagnetic

Speed: 0–6 km/h

Seat dimensions: Width 50 cm × Depth 47 cm × Thickness 4.5 cm

Seat height from the ground: 52 cm
Front wheels: 8" honeycomb rubber wheels
Rear wheels: 12" honeycomb rubber wheels
Armrest height from the ground: 73 cm
Internal width between armrests: 56.5 cm
Braking distance: ≤ 1000 mm
Minimum turning radius: ≤ 900 mm
Driving range: Up to 20 km
Static stability: $\geq 6^\circ$
Dynamic stability: $\geq 6^\circ$
Maximum slope angle: 6°
Climbing ability: $\geq 6^\circ$
Maximum user weight: 220 kg
Net weight (without battery): 22 kg
Gross weight: 33 kg
Additional Information
Can the wheelchair be used as a seat in a motor vehicle? No
NOTE: Product dimensions may vary by ± 2 cm.

This symbol indicates the maximum user weight.

ASSEMBLY

Remove the wheelchair and all remaining components from the packaging.


Install the controller on either the right or left armrest.

Hold the push handle with one hand and the seat frame with the other, then unfold the wheelchair by pulling the frame outward.

Lock the clamp by placing the U-shaped clamp onto the latch and pressing the lever downward.

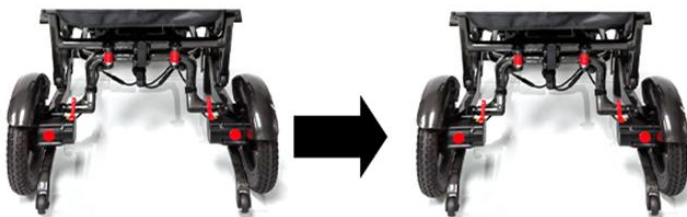Lock the armrest by pressing it down and then rotating the small lever located on the side of the armrest upward. (Depending on the model, the release mechanism may be one of two different designs.)

CONTROLLER

- 1.LED Display
- 2.ON/OFF Button
- 3.Speed Decrease Button
- 4.Horn
- 5.Battery Level Indicator
- 6.Speed Level Indicator
7. Speed Increase Button

Power On

Press the ON/OFF button. The battery level indicators will illuminate sequentially from left to right.

Power Off

Press the ON/OFF button. All LED indicators will turn off.

Speed Adjustment

The speed of the electric wheelchair can be adjusted according to the user's preferences and the surrounding environment. The speed level indicator shows the currently selected maximum speed. Use the Speed Increase and Speed Decrease buttons to adjust the speed.

Horn

Press the Horn button to activate the audible warning signal. The horn will sound for as long as the button is pressed.

Battery Level Indicator

When the electric wheelchair is switched on, the battery level indicator displays the approximate remaining battery capacity. When all LED indicators are illuminated, the battery is fully charged.

SAFETY INSTRUCTIONS

Do not carry passengers.

Do not use the wheelchair on steep slopes.

Do not operate the wheelchair under the influence of alcohol.

Do not tow other vehicles.

Do not use or operate handheld communication devices, such as two-way radios or mobile phones, while driving.

Do not use the wheelchair without first reading the user manual.

Do not exceed the maximum safe incline angle of 13°.

Do not make sudden joystick movements while driving on slopes.

Do not carry passengers or exceed the maximum user weight while driving.

Do not press the controller ON/OFF button while driving.

Do not drive on soft terrain such as loose soil, gravel, or tall grass.

Do not attempt to climb steps higher than 50 mm unless a suitable ramp is provided.

Do not disassemble the wheelchair unless the electromagnetic brakes are engaged and the controller is switched off.

Do not operate the wheelchair when no one is seated in it.

Do not drive the wheelchair on the roadway. Use only sidewalks and designated pedestrian routes.

Do not change direction or reverse until the wheelchair has come to a complete stop.

Always select the lowest speed when descending a slope.

Always fasten the seat belt.

Always keep your feet on the footrests while driving.

Do not stand on the footrests.

Before driving, make sure the batteries are fully charged.

Charge the wheelchair in a well-ventilated area.

If the center wheel bolt becomes loose, replace it with a standard bolt using an appropriate wrench.

When descending a slope, move the seat as far forward as possible.

The wheelchair is designed to transport one person only.

WARNINGS AND RECOMMENDATIONS

Do not use the electric wheelchair if the motor drive is disengaged.

Do not disassemble the motor, as this may affect the wheelchair's performance and cause damage.

Precautions When Driving on Slopes

Do not climb slopes steeper than 12°.

Do not reverse while driving uphill, as this may cause the wheelchair to tip over.

Do not drive across side slopes exceeding 3°, as this significantly increases the risk of tipping over.

Note: The controller can be mounted on either the left or right armrest according to the user's needs. The controller mounting tube is secured with a knob located on the controller bracket and can be adjusted forward and backward as required.


JOYSTICK AND CONTROLLER OPERATING INSTRUCTIONS

Installation

Attach the controller to the controller mounting bracket.

Adjust the length of the controller mounting tube, then tighten the nut securely.

Connect the cables to the left motor, right motor, and battery.

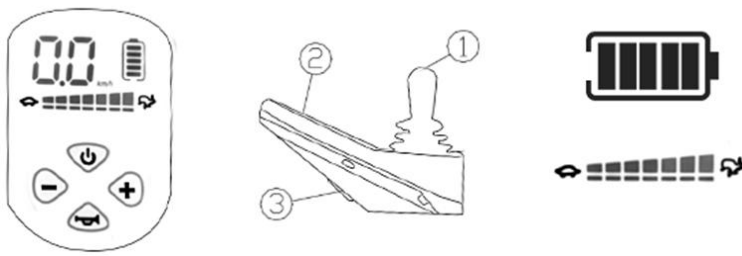CONTROLLER OPERATION

Avoid striking the controller, especially the joystick. Do not hit or force the joystick while driving. Before operating the wheelchair, make sure that the control system is securely installed and functioning properly. Do not store the controller or its components in a damp environment for extended periods. If the controller comes into contact with food or beverages, clean it immediately.

Cleaning

Clean the controller and joystick using a soft cloth dampened with a diluted disinfectant solution. Take extra care when cleaning the joystick. Do not use abrasive cleaners or other aggressive cleaning agents.

REAR WHEEL CLUTCH LEVER

Electric Mode

When the wheelchair stops in electric mode, the electromagnetic braking system is automatically engaged, preventing the wheelchair from being pushed manually. To enable manual movement, move the red clutch lever to Manual Mode. In this mode, the electric drive is disengaged.

Manual Mode

To switch the wheelchair to manual mode, simply move the red clutch lever to the appropriate position.

Note: It is recommended to keep the wheelchair in Electric Mode whenever possible. The clutch function is intended only to allow manual pushing of the wheelchair when necessary, for example in shops or crowded areas.

Manually Pushing the Wheelchair

Move the motor clutch levers on both the left and right sides to the Manual Mode position. The wheelchair can then be pushed manually.

FOLDING THE WHEELCHAIR

Pull the small lever located on the outside of each armrest downward to unlock the armrests. (On some models, the armrests are released by pressing buttons instead of using levers.)

Release the rear locking mechanisms and lower the U-shaped locking bracket.

Hold the push handle with one hand and the seat frame with the other, then fold the wheelchair by bringing both parts together.

Armrests

Pull the small lever located on the outside of each armrest downward to unlock the armrests. Once unlocked, the armrests can be tilted backward.

For users who are unable to move their legs or operate the armrests independently, assistance from a caregiver is recommended.

Important: Anyone assisting the wheelchair user should be properly trained in the safe operation of the wheelchair.

BATTERY

Removing the Battery

Disconnect the power cable. Press the locking latch to release the battery from the frame, then slide the battery out of its compartment.

Battery Use

The electric wheelchair is equipped with a 24 V sealed, maintenance-free lithium-ion battery designed for deep discharge cycles.

Although it may resemble a car battery, it must never be replaced with one, as automotive batteries are not designed for deep-cycle operation and are not safe for use in electric wheelchairs.

BATTERY CHARGING AND DISCHARGING

Before first use, fully charge the battery to achieve approximately 90% of its rated performance.

During the first few uses around your home or yard, drive slowly and avoid traveling too far until you become familiar with the wheelchair's controls. Then fully discharge the battery.

After the battery has been completely discharged, charge it for 8–14 hours before using the wheelchair again. This procedure allows the battery to reach more than 90% of its full performance.

After 4–5 complete charge and discharge cycles, the battery will reach its optimum performance and service life.

Important Battery Information

Regular charging ensures reliable operation and extends battery life. Whenever possible, keep the battery fully charged.

Infrequent or partial charging may permanently damage the battery, reduce its capacity, and shorten its service life.

If the wheelchair is not used regularly, charge the battery at least once a week.

CHARGER

The charger is one of the most important components of the electric wheelchair and provides safe, fast, and convenient battery charging.

Maintenance Recommendations

The motor carbon brushes should be inspected every 6 months by an authorized service center.

Protect the controller and joystick from moisture, as it may cause damage.

Regularly inspect the battery cables for signs of corrosion.

Clean the wheelchair using a damp cloth. Do not pour water over the wheelchair or wash it under running water.

ELECTRIC WHEELCHAIR TROUBLESHOOTING AND SELF-CHECK

If a malfunction occurs while operating the electric wheelchair, switch the power OFF before performing any inspection.

A complete power failure is indicated when all LED indicators on the control panel are off.

Inspection

Check that the controller plug is properly connected and not loose.

Check that the cable connecting the controller to the battery compartment is securely connected. If necessary, disconnect and reconnect the connector by holding the connector itself rather than the cable to avoid damaging the wires.

If the power is not restored after completing the above checks, or if you have any questions, contact an authorized service center immediately.

Note: The controller is equipped with a diagnostic system that continuously monitors the operation of both the controller and the motors.

If a fault is detected, the corresponding error code will be displayed on the controller screen.

IMPORTANT NOTICE – Serious incidents

In the event of a “serious incident” related to the medical device that directly or indirectly led, could have led, or may lead to any of the following events:

a) death of a patient, user, or other person, or

b) temporary or permanent deterioration of the health of a patient, user, or other person, or

c) a serious public health threat,

the above “serious incident” must be reported to the manufacturer and to the competent authority of the Member State in which the user or patient resides.

NOTE: If pain, allergic reactions, or other disturbing or unclear symptoms related to the use of the medical device occur, consult a healthcare professional.

NOTE: It is prohibited to use the product in any manner other than its intended purpose.

CLEANING AND MAINTENANCE

User-performed maintenance includes cleaning, checking the tires, the battery, and the charging process. More advanced maintenance should be carried out by service personnel. It is recommended to inspect the wheelchair every two months and perform a full technical inspection once a year.

Cleaning the wheelchair

Regularly remove dust from the wheelchair surfaces. Use wax polish or a clean, soft cloth for cleaning.

Tires

Regularly check the wear of pneumatic tires and replace them when the tread depth reaches 1 mm.

Battery

The battery should be charged regularly. Fully discharging the battery before recharging is not recommended, as it shortens battery life.

STERILIZATION AND DISINFECTION

The wheelchair is not a sterile medical device. During daily use, it should be kept clean.

If the wheelchair becomes contaminated with blood, secretions, excretions, or other substances, wipe it with a disinfectant based on ammonium salts at a concentration of 1000–2000 mg/L, leave it for 30 minutes, then wipe with a clean, dry cloth and allow it to dry.

WHEELCHAIR DAMAGE AND SELF-CHECK

If any abnormalities occur during use, switch off the power and check the wheelchair.

Symptom: complete loss of power and the controller indicators are off.

Step 1: Check whether the controller plug is loose or disconnected.

Step 2: Check the connection between the battery box and the controller. If necessary, reconnect the plug by holding the connector, not the cable, to avoid damage to the wiring.

If power is not restored after performing the above steps or if other concerns arise, contact service immediately.

STORAGE AND TRANSPORT CONDITIONS

The electric wheelchair should be stored and transported under the following conditions:

Temperature range: 5°C – 40°C Relative humidity: ≤ 80% Atmospheric pressure: 860 hPa – 1060 hPa

Rated voltage Charger operating voltage: AC 220 V, 50 Hz Internal power supply voltage: DC 24 V

DISPOSAL: At the end of its service life, the product must be disposed of in accordance with applicable local, national and EU regulations concerning waste electrical and electronic equipment (WEEE Directive). The device must not be disposed of together with unsorted municipal waste. The product contains recyclable components, such as plastics, electronic components and batteries. Particular attention should be paid to the battery, which contains substances requiring appropriate handling and should be handed over to specialised collection points. For information on proper disposal methods, contact local authorities or an authorised seller. **NOTE:** Incorrect disposal may pose a risk to the environment and human health.

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Model:
Serial number:

Point of sale stamp

and a legible signature of the Seller:

1. ANTAR MEDICAL Sp z o.o, 03-068 Warsaw, ul. Zawisłańska 43, hereby grants a 12-month warranty on the product, starting from the date of purchase by the Buyer.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, please contact your retailer/distributor or (if it is impossible) send it to ANTAR (as specified in point 1).
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum allowed weight of the user, as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts which are subject to normal wear or tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic hazard.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.
10. The warranty does not cover parts and components subject to natural wear and tear, such as wheels, tires, rubber pads, fuses, etc., the wear of which results from normal use of the product. The warranty period for batteries is 6 months.

ANTAR MEDICAL Sp. z o.o
03-068 Warsaw, ul. Zawisłańska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net

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