AT52338

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EN

Thank you for purchasing our product. Please read this user manual carefully before use. Do not operate the product before re ading and fully understanding these instructions. If you do not understand any warnings, cautions, or recommendations, contact a health care professional or your dealer to avoid product damage or personal injury. WARNING: All parts must be checked for damage. Before first use, thoroughly inspect the technical condition of the product and ensure that all parts are complete and properly assembled (if a pplicable, pay particular attention to screws, nuts, rivets and other fastening and structural components). Make sure that the product is not defective or damaged. If any damage or missing parts are found, DO NOT use the product. For further information, contact the seller.

TARGET USER GROUP

This product is intended for individuals with conditions, dysfunctions, or injuries requiring treatment, rehabilitation, or mobility support (see the intended use section of this manual).

The product may be purchased independently by the user or based on the recommendation of a physician, therapist, or other spe cialist. In both cases, the available sizes, required features, product variants, indications and contraindications, and the manufacturer's information must be taken into account.

INTENDED USE: The electric wheelchair is designed as a mobility vehicle for individuals with disabilities. It is powered by an electric motor, with a design that limits the speed to a safe walking pace.

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post-injury or post-surgical conditions requiring reduced load on the lower limbs,

neurological conditions affecting walking ability,

elderly individuals with limited mobility.

The wheelchair is intended for use by a single person who is able to operate the device independently or under caregiver supe rvision.

CONTRAINDICATIONS

Physical or mental limitations (e.g., impaired vision) that prevent safe operation of the product.

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GENERAL SAFETY INSTRUCTIONS

MANDATORY! Read this manual carefully before first use and ensure you understand all instructions.

The electric wheelchair is a medical device designed to improve user mobility. Safe operation requires caution, attention, an d adjustment of driving behavior to environmental conditions.

Drive slowly until you are familiar with all functions and features.

Reduce speed in crowded areas.

Pay special attention to pedestrians and obstacles.

Use the wheelchair only in accordance with local regulations.

PROHIBITED! Keep hands away from tires and wheels while driving. Loose clothing may be caught in moving parts.

USAGE WARNINGS

WARNING! Do not exceed the maximum load of 130 kg (user + cargo).

WARNING! The wheelchair is intended for one person only. Do not carry passengers, including children.

WARNING! Do not use the wheelchair on stairs or escalators.

WARNING! Do not attempt to overcome obstacles higher than specified and always approach them head-on.

WARNING! Do not use the wheelchair in heavy rain, high humidity, or extreme temperatures.

WARNING! Surfaces of the wheelchair may become hot – use caution.

DRIVING SAFETY

Always adjust your speed to road conditions. Pay attention to pedestrians and adapt your driving style accordingly. In crowde d areas, always move at reduced speed.

Slow down before turning – sharp turns at high speed may cause tipping.

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WARNING! Do not lean excessively forward or sideways – this may cause loss of balance and tipping.

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SAFETY ON SLOPES AND RAMPS

Do not exceed the maximum allowable slope.

Do not drive across slopes or diagonally.

Use minimum speed when descending.

Maintain steady speed when ascending.

When driving in crowded indoor environments (e.g., supermarkets), maintain low speed to control braking distance. When turnin g, do not exceed 2 km/h.

Do not use the wheelchair in rain and avoid driving through water. If unavoidable, water depth must not exceed 20 mm and spee d must not exceed 1 km/h.

When not in use, switch off the power using the key and remove it to prevent unintended operation by children or unauthorized persons.

PRE-USE CHECK

Before each use, check: battery condition, brake operation, electrical connections, tire condition, proper operation of all functions. If any irregularities are detected, do not use the wheelchair and contact service.

BATTERY AND ELECTRICAL SAFETY

WARNING! Do not allow the battery to come into contact with metal objects – risk of short circuit and electric shock.
WARNING! Use only recommended batteries and chargers.
WARNING! Do not charge the battery near flammable materials.
WARNING! Do not allow children near the wheelchair during charging.

CHILDREN AND BYSTANDERS

WARNING! Do not allow children to operate the product.
WARNING! Do not leave the wheelchair unattended near children.
WARNING! Do not carry passengers – this may lead to instability.

MODIFICATIONS AND PARTS

WARNING! Do not modify the wheelchair or use non-original parts or accessories.
WARNING! Do not lift the wheelchair by removable parts (e.g., armrests, seat).

EMERGENCY SITUATIONS

WARNING! If the wheelchair does not stop properly, immediately use the brake, switch off power, and contact service.
WARNING! After an accident or damage, the wheelchair must be inspected by qualified service before reuse.

ELECTROMAGNETIC INTERFERENCE (EMI/RFI)

WARNING! Devices such as mobile phones or radios may cause unintended movement.
WARNING! Do not use radio devices while driving and avoid strong signal sources.

TRANSPORT IN A MOTOR VEHICLE

During transport, the electric wheelchair must be carried without the user, in the trunk or cargo area, with batteries removed and properly secured against movement.
Removable parts should be detached or firmly secured.
WARNING! Do not sit on the wheelchair during transport.
WARNING! The wheelchair and batteries must be properly secured.
WARNING! Do not transport batteries with flammable materials.

COMPONENT DESCRIPTION

1.Ignition switch (key)
2.Throttle lever
3.Battery display
4.Brake lever
5.Forward/reverse switch

6.Horn button
7.Waterproof wiring harness
8.Folding steering column
9.Headlights / horn
10.Front fender
11.Front wheel
12.Foot platform
13.Lithium battery
14.Aluminum alloy frame
15.Seat cushion
16.Backrest
17.Lock latch
18.Basket
19.Rear fender
20.Reflector
21.Drum brake
22.Anti-tip wheel
23.Motor

FOLDING / UNFOLDING

1.Open the red latch to fold/unfold the seat
2.Open the black safety lock to fold/unfold the steering column
3.Remove the lithium battery before folding
4.The folding column can be adjusted up and down

ACCESSORIES DESCRIPTION

Waterproof cable connector

The waterproof plug and connector ensure safe and reliable electrical and signal connections. They are highly resistant to moisture and easy to replace.

Electric control switch

The electric switch controls lights, horn, and forward/reverse movement. The buttons are designed to be anti-slip, symmetrical, and easy to operate.

Red button: forward/reverse switch

Green button: light switch

Brake lever (black)
The brake lever is made of aluminum alloy – lightweight, corrosion-resistant, and durable. In combination with the braking system, it ensures quick and effective braking.

Throttle control
Provides comfortable driving and precise control of the wheelchair. The ergonomic grip is adapted to the hand and features an anti-slip surface. It displays driving data, allowing the user to monitor operating parameters in real time.

Ignition switch (key)
Position 1: power off
Position 2: power on
Position 3: power on + lights

PACKAGE CONTENTS: Wheelchair, battery, charger, cushion, basket, user manual, toolkit
TECHNICAL SPECIFICATIONS
Unfolded dimensions (length × width × height): 1010 × 530 × 790 mm
Folded dimensions (length × width × height): 1010 × 530 × 370 mm
Seat dimensions (W × L × T): 420 × 320 × 30 mm
Seat height from ground: 460 mm
Backrest height: 350 mm
Motor: 36V 300W (2 × brushless)
Battery: 36V 6Ah lithium
Charger: 42V 2A
Controller: 36V 15A (brushless)
Handlebar: height adjustable
Speed: ≤ 6 km/h
Range: ≤ 18 km
Brake: electronic + mechanical
Braking distance: 1000 mm
Minimum turning radius: ≤ 1200 mm
Ground clearance: 95 mm
Obstacle climbing ability: 30 mm
Maximum slope: ≤ 15°
Static stability: ≥ 9°, < 15°
Dynamic stability: ≥ 6°, < 10°
Lighting: front light + rear reflector
Basket: front and rear
Front basket dimensions: 340 × 260 × 250 mm
Front basket weight: 0.8 kg
Front basket max load: 25 kg
Rear basket dimensions: 340 × 295 × 150 mm
Rear basket weight: 1.15 kg
Rear basket max load: 25 kg
Front wheel: 8" solid PU tires, nylon hub
Rear wheel: 8" solid PU tires, aluminum alloy hub
Maximum load: ≤ 130 kg
Net weight: 21.8 kg
Gross weight: 31 kg
MATERIAL COMPOSITION
Frame: aluminum alloy
Motor: metal
Battery: lithium battery
Battery housing: plastic
Charger: plastic

Armrests: none
Seat: PVC cover + PU foam filling
Backrest: PVC cover + PU foam filling
Fenders: ABS
Front basket: nylon fabric
Rear basket: steel
Handle grip: PU foam
Plastic parts: ABS

NOTE: Individual product dimensions may vary by ±2 cm

This symbol indicates the maximum user weight.

FRONT BASKET INSTALLATION

Place the wheelchair on a stable, level surface. Turn off the power and make sure that the wheelchair cannot move by itself.
Prepare the front basket, mounting bracket, fixing screws and the appropriate Allen key.
Place the mounting bracket at the front of the steering column, below the handlebar. The bracket should be positioned centrally on the front side of the column.

When positioning the bracket, make sure that the cables located near the handlebar are not pinched, bent or pulled tight.
Insert the fixing screws into the holes of the mounting bracket. Tighten them gradually using the Allen key, alternately on both sides, so that the bracket fits evenly against the steering column.
After initially tightening the screws, check the position of the bracket. It should be mounted straight and securely, and it should not touch or block the control cables.
Tighten all fixing screws until the bracket is securely attached. Do not overtighten the screws to avoid damaging the mounting parts.

Hold the front basket and position its rear part opposite the installed bracket.
Slide the basket onto the mounting bracket from above, so that the basket hooks fit into the corresponding places on the bracket.

Press the basket down until it is properly seated on the bracket.

After installation, gently move the basket to check whether it is securely attached and does not come off the bracket.
Before driving, make sure that the basket does not limit the steering movement, block the cables or affect the operation of the control elements.

BATTERY INSTRUCTIONS

Improper use will shorten the battery's service life and may lead to overheating, fire, or even explosion. Please read the following instructions carefully:

Do not disassemble or dismantle the battery. It contains protective circuit components – damage to these parts may lead to malfunction.

It is strictly prohibited to short-circuit the positive and negative terminals. A short circuit may cause overheating, fire, or other hazards.

Do not burn or heat the battery. This may melt the insulation and result in fire, explosion, or other dangers.

Do not immerse the battery in water or use it in a humid environment. Water may cause a short circuit in the protection system, leading to overheating or other risks.

Charge the battery only with a compatible charger. Using an incorrect charger may cause abnormal chemical reactions, overheating, or fire.

Do not strike the battery with a hammer or other hard objects, and do not cut it with sharp tools. A damaged battery may cause short circuits and fire.

If the battery is damaged or leaking and the substance comes into contact with the eyes, do not rub them. Rinse immediately with clean water and seek medical attention.

Do not mix lithium-ion batteries from different manufacturers or with different capacities, as this may lead to overheating, fire, or other hazards.

After use, fully charge the battery within 3 days. If not used for a longer period, recharge it every two months to prevent damage.

Do not exceed the recommended charging time (usually up to 10 hours) to avoid overheating. If the battery is not fully charged after this time, stop charging.

Keep the battery out of reach of children.

Do not store the battery in a microwave or in a high-pressure environment, as this may result in serious danger.

TROUBLESHOOTING

Speed controller malfunction or significant reduction in maximum speed

Causes:

Loose battery cable.

Low battery charge.

Damaged controller spring.

Solutions:

Tighten the battery cable

Fully recharge the battery.

Replace the spring or controller.

Charger does not charge property.

Causes:

Damaged power cable or loose socket

Blown charger fuse.

Disconnected battery terminals.

Solutions:

Repair the cable or secure the connection

Replace the fuse with the same type.

Reconnect all cables.

Charging is normal but range is insufficient

Causes:

Low tire pressure or overly tight brake adjustment.

Prolonged high-speed driving, uphill travel, headwinds, frequent kraking.

Worn or damaged battery.

Solutions:

Inflate tires and properly adjust the brake

Avoid long high-speed driving and frequent braking (increased power consumption uphill and against wind is normal).

Replace or repair the battery.

Drive motor does not function property.

Causes:

Loose or disconnected cables (battery, motor, controller).
Damaged or contaminated contro ller.
Solutions:
Check and reconnect all cables.
Keep the controller clean and free from contamination.

IMPORTANT NOTICE

In the event of a serious incident related to the device that directly or indirectly led, could have led, or may lead to:
a) death of a patient, user, or other person, or
b) temporary or permanent deterioration of health of a patient, user, or other person, or
c) a serious public health threat
such incident must be reported to the manufacturer and the competent authority of the Member State where the user or patient resides.
NOTE: If pain, allergic reactions, or other concerning or unclear symptoms occur during use of the medical device, consult a healthcare professional.
NOTE: The product must not be used in any manner other than its intended purpose.

CLEANING AND MAINTENANCE

Cleaning
Clean the product as needed using a soft, damp cloth and mild detergents. After cleaning, wipe with clean water and dry thoro ughly.
Do not use pressure washers or large amounts of water, as this may damage electrical components. Cleaning should be performed at temperatures above 0°C.
Maintenance
Although the scooter requires minimal servicing, regular checks are necessary:
Daily: check battery charge level, horn operation, and cleanliness of the battery compartment
Weekly: check throttle lever, brake operation, and stability of wheels and connections
Monthly: check tire wear and proper driving and braking performance
Periodically: inspect wiring, connections, and all structural components
All repairs, adjustments, and battery replacement must be performed only by qualified service personnel.

STORAGE

Store the product in a dry place. Do not expose it to direct sunlight.

DISPOSAL

At the end of its service life, the scooter must be disposed of in accordance with applicable local, national, and EU regulat ions regarding waste electrical and electronic equipment (WEEE Directive).
Do not dispose of the device with unsorted municipal waste. The scooter contains recyclable materials such as plastics, electronic components, and batteries.
Special attention must be paid to the battery, which contains substances requiring proper handling (including lithium compoun ds) and must be disposed of at designated collection points.
For information on proper disposal methods, contact local authorities or an authorized dealer.
WARNING: Improper disposal may pose a risk to the environment and human health.

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Serial number:

Point of sale stamp

and a legible signature of the Seller:

1. ANTAR MEDICAL Sp. z o.o. ,03-068 Warsaw, ul. Zawiślańska 43, hereby grants a 12-month warranty on the product, starting from the date of purchase by the Buyer.
2. During the warranty period ANTAR undertakes to carry out the repairs necessary to restore proper functioning of the product, within 14 days from the date of receipt of the product by the service centre, together with the proof of purchase and the warranty card.
3. ANTAR reserves the right to replace the product if the repair is deemed uneconomic. Repair or replacement of faulty products does not extend the warranty period.
4. If the product is found to be malfunctioning, please contact your retailer/distributor or (if it is impossible) send it to ANTAR (as specified in point 1).
5. Before use, read the instructions carefully and follow them. If the product is used improperly and contrary to recommendations, the warranty will be irrevocably lost.
6. The warranty does not cover damage caused by external mechanical force, load exceeding the maximum allowed weight of the user, as described in the manual, as well as contact of the product with alcohol, greasy substances or petrol. The warranty does not cover parts which are subject to normal wear or tear during use, such as the bearing, rubber parts and others.
7. The product delivered to the service in connection with the complaint cannot cause an epidemic hazard.
8. THIS WARRANTY IS ONLY VALID TOGETHER WITH THE PROOF OF PURCHASE (BILL, RECEIPT, INVOICE). IN CASE OF A COMPLAINT, A PROOF OF PURCHASE MUST BE ATTACHED.
9. Complaints submitted without the proof of purchase and a warranty card with the product serial number will not be considered.
10. The warranty does not cover parts and components subject to natural wear and tear, such as wheels, tires, rubber pads, fuses, etc., the wear of which results from normal use of the product. The warranty period for batteries is 6 months.

ANTAR MEDICAL Sp. z o.o.
03-068 Warsaw, ul. Zawiślańska 43
Tel. 22 518 36 00, Fax 22 518 36 30
www.antar.net

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En

Thank you for purchasing our product. Please read the user manual carefully. Do not use the product before reading and understanding these instructions. If you do not understand the warnings, notes, and recommendations, contact a healthcare professional or the seller to avoid product damage or personal injury. WARNING: All parts must be checked for damage. Before first use, thoroughly inspect the technical condition of the product and ensure that all parts are complete and properly assembled (if applicable, pay particular attention to screws, nuts, rivets and other fastening and structural components). Make sure that the product is not defective or damaged. If any damage or missing parts are found, DO NOT use the product. For further information, contact the seller.

TARGET PATIENT GROUP: Individuals suffering from conditions, dysfunctions, or injuries for which this product is intended for treatment/rehabilitation or compensation (see the section on intended use in this manual). The product may be purchased independently by the user or based on the recommendation of a physician, therapist, or other specialist. In both cases, the available sizes/required functions and product variants, indications and contraindications for use, as well as the information provided by the manufacturer must be taken into account.

INTENDED USE

An electric wheelchair is a vehicle designed for the mobility of a person with disabilities. It is powered by a motor, the design of which limits the driving speed to walking pace.

The AT52339 mobility scooter has been designed for users with a maximum weight of up to 120 kg, who move in a seated position. Safe use of the device requires an appropriate level of motor skills, vision, and cognitive ability; therefore, consultation with a physician before use is recommended.

The AT52339 model is intended for both indoor and outdoor use, within a wide temperature range from -20°C to +40°C. It can be used during light rain; however, prolonged use in heavy rain or snow is not recommended.

Thanks to relatively large wheels, the scooter allows driving on unpaved surfaces, moderately uneven terrain, and overcoming obstacles up to 4 cm in height. The maximum safe incline is 9° (16%) and must not be exceeded.

The vehicle is easy to operate, and the braking system ensures smooth and quick stopping after releasing the acceleration lever. Under standard conditions, the scooter range is approximately 15 km; however, this may decrease when driving uphill, on uneven terrain, or when frequently overcoming obstacles.

Before first use, carefully read the user manual and perform a visual inspection of the device to ensure that it has no visible damage. In case of any doubts, contact the seller.

GENERAL SAFETY RULES

MANDATORY! Before first use, carefully read this manual and make sure you understand all the rules for using the AT52339 mobility scooter.

The AT52339 mobility scooter is a medical device intended to improve user mobility. Its safe use requires caution, concentration, and adaptation of driving style to environmental conditions.

Drive slowly until you become familiar with all functions and features of the scooter.

In crowded areas, move at reduced speed.

Exercise particular caution with pedestrians and obstacles.

Use the scooter only in accordance with local regulations.

PROHIBITED! While driving, keep your hands away from tires and wheels. Remember that loose clothing may be caught by the tires and wheels.

WARNINGS REGARDING USE

WARNING! Do not exceed the maximum permissible load of 120 kg (user + luggage).

WARNING! The scooter is intended for one person only. Do not carry passengers, including children.

WARNING! Do not use the scooter on stairs or escalators.

WARNING! Do not overcome obstacles higher than specified and always approach them frontally.

WARNING! Do not use the scooter in heavy rain, high humidity, or extreme temperatures.

WARNING! Scooter surfaces may become hot – exercise caution.

DRIVING SAFETY

Always adjust your speed to road conditions. Pay particular attention to pedestrians and adapt your driving accordingly. In crowded places, always move at half speed.

Reduce speed before turning – turning too quickly may cause the vehicle to tip over.

Maintain a stable body position while driving.

Avoid sudden maneuvers and abrupt braking.

WARNING! Do not lean excessively forward or sideways – this may result in loss of balance and tipping of the scooter.

PROHIBITED! Do not place hands or feet near moving parts (wheels, mechanisms).

SAFETY ON SLOPES AND RAMPS

Do not exceed the permissible incline.

Do not travel across or diagonally on slopes.

When descending, set the minimum speed.

When ascending, maintain a constant speed.

WARNING! Never use freewheel mode on a slope.

PRE-USE INSPECTION

Before each ride, check: battery condition, brake operation, electrical connections, tire condition, and proper functioning of all features. If any irregularities are detected, do not use the scooter and contact service.

BATTERIES AND ELECTRICAL SAFETY

WARNING! Do not allow the battery to come into contact with metal objects – risk of short circuit and electric shock.

WARNING! Use only recommended batteries and chargers.

WARNING! Do not charge the battery near flammable materials.

WARNING! Do not allow children near the scooter during charging.

CHILDREN AND BYSTANDERS

WARNING! Do not allow children to operate the AT52339 mobility scooter.

WARNING! Do not leave the scooter unattended in the presence of children.

WARNING! Do not carry passengers. This may cause loss of stability.

MODIFICATIONS AND PARTS

WARNING! Do not modify the scooter or use non-original parts and accessories.

WARNING! Do not lift the scooter by detachable components (e.g., armrests, seat).

EMERGENCY SITUATIONS

WARNING! If the scooter does not stop properly, immediately use the brake, switch off the power, and contact service.

WARNING! After an accident or damage, the scooter must be checked by qualified service before reuse.

ELECTROMAGNETIC INTERFERENCE (EMI/RFI)

WARNING! Devices such as mobile phones or radios may cause unintended scooter movement.

WARNING! Do not use radio devices while driving and avoid staying near strong signal sources.

TRANSPORT IN A MOTOR VEHICLE

During transport, the AT52339 electric mobility scooter must be transported without the user, in the trunk or cargo space of the vehicle, with batteries removed and properly secured against movement.

Detachable elements should be removed or firmly secured.

WARNING! Do not sit on the scooter during transport.

WARNING! The scooter and batteries must be properly secured.

WARNING! Do not transport batteries with flammable materials.

INDICATIONS

The product is intended for use in cases of:

difficulty walking independently or rapid fatigue during movement,

musculoskeletal disorders (e.g., joint diseases, injuries, degeneration),

reduced muscle strength or balance disorders,

conditions after injuries or surgeries requiring relief of the lower limbs,

neurological diseases affecting walking ability,

advanced age associated with limited mobility.

The scooter is intended for use by one person who is able to operate the device independently or use it under caregiver super vision.

CONTRAINDICATIONS: Physical or mental limitations (e.g., visual impairment) that prevent safe handling of the product.

Lack of ability to safely operate the vehicle in outdoor environments (e.g. lack of orientation in pedestrian or road traffic).

Uncontrolled body movements (e.g. tremors, spasticity) that may hinder operation.

Episodes of loss of consciousness, dizziness, or other conditions that may occur suddenly أثناء driving.

Inability to maintain a stable seated position.

Significant limitations in reaction time and motor coordination.

Use without supervision in the case of persons requiring constant care.

LIST OF COMPONENTS: Electric wheelchair, battery, charger, cushion, basket, user manual, tool kit.

TECHNICAL SPECIFICATION

Dimensions when unfolded (length x width x height): 970 × 530 × 900 mm

Dimensions when folded (length x width x height): 970 × 530 × 290 mm

Seat dimensions (W × L × T): 390 × 340 × 30 mm

Motor: 24V 250W (2 × brushless)

Battery: 24V 8A lithium, removable

Battery weight: 2 kg

Charger: 24V 2A

Controller: 24V 60A brushless

Steering column: height adjustment

Speed: adjustable 0–6 km/h

Brake: electromagnetic ABS

Braking distance: ≤ 1000 mm

Minimum turning radius: ≤ 1500 mm

Ground clearance: 80 mm

Obstacle climbing ability: 40 mm

Backrest height: 830 mm

Armrest height: 720 mm

Distance between armrests: 435 mm

Static stability: ≥ 9°, < 15°

Dynamic stability: ≥ 6°, < 10°

Lighting: front and rear light

Basket: front basket

Basket dimensions: 310 × 230 × 210 mm

Maximum basket load: 25 kg

Range: 15 km

Front wheel: 9" pneumatic tires, aluminum alloy hub

Rear wheel: 9" pneumatic tires, aluminum alloy hub

Maximum load: 120 kg

Net weight: 18.6 kg
Gross weight: 28 kg
MATERIAL COMPOSITION
Frame: aluminum alloy
Motor: metal
Battery: lithium battery
Battery housing: plastic
Charger: plastic
Armrests: PU foam
Seat: PU foam
Backrest: PU foam
Mudguards: plastic
Front basket: nylon fabric
Rear basket: none
Steering system: ABS
Plastic components: ABS
NOTE! INDIVIDUAL PRODUCT DIMENSIONS MAY VARY WITHIN A TOLERANCE OF +-2 CM

This symbol indicates the maximum user weight

1.Control panel
The control panel includes the display, speed adjustment knobs, key, lights, and horn.
2.Charging port
For charging.
3.USB socket
For mobile phone, etc., output voltage 5 V.
4.Hook
Used for hanging items, maximum load 3 kg.
5.Steering column adjustment knob
Used to adjust the angle of the steering column.
6.Front lights
7.Removable battery
The battery can be removed.
8.Basket / storage compartment
Maximum load 6 kg.
9.Armrest adjustment knob
After loosening the screws, the armrest can be adjusted outward. It should then be tightened again.
10.Rear armrest folding button
Press the button to fold the armrest backward.
11.Seat cushion
12.Armrest
The seat width can be adjusted, and the backrest can be reclined backward.
13.Backrest
Can be folded forward.

14.Folding handle
Special care must be taken when using the folding handle. The seat folding mechanism is locked by the handle. Before use, pre ss the seat with your hand to make sure that it is properly locked. The condition of the folding handle should be checked regularly. More information can be found below.
15.Rear light
16.Rear bumper
17.Anti-tip wheel
18.Mudguard


1.Power switch
When the power is switched off, the scooter cannot move electrically.
Turn the key clockwise to switch on the power and enable driving.
2.Battery display
After switching on the power, the battery level indicator shows its status.
The green area indicates a high charge level, while the red area indicates a low level and the need for immediate charging.
3.Speed control
The speed control allows the maximum forward driving speed to be set in the range from 0 to 6 km/h.
Turning the knob fully to the left sets the lowest speed, and fully to the right sets the highest speed.
4.Front and rear light switch
Turn on the lights when driving at night or in conditions of limited visibility.
5.Horn button
Pressing the button emits an audible signal.
To prevent accidents and ensure safe driving, use it when necessary.
6.Forward and reverse drive lever
Push the right lever forward to drive forward; push the left lever forward to drive backward.
When the lever is released, it automatically returns to the “middle” position, which activates the electromagnetic brake and stops the scooter.
7.Steering handlebar
While driving, hold both hands firmly on the handlebar.

OPERATING PRINCIPLE OF THE FOLDING LEVER

1.Folding handle
2.Spring
3. Pin/bolt

BRAKE CONTROL LEVER

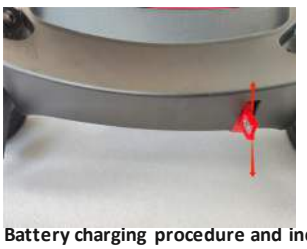Battery charging procedure and indication

The charger should be connected to the charging port located on the battery housing. The charger has protection against incor rect connection (it cannot be inserted incorrectly).
During charging, a red indicator light is on, indicating the charging process. Once the battery is fully charged, a green ind icator light comes on, indicating the end of charging.
WARNING! Never connect the charger incorrectly, as this may damage the charging port.

Overload protection and fuse


The battery housing contains overload protection. This is an important safety feature that protects the battery from damage b y controlling excessive current flow during charging.
The fuse is integrated with the overload protection. In the event of a short circuit, the fuse blows, protecting the entire e lectrical system.

MANUAL AND AUTOMATIC SWITCHING OF THE BRAKE LEVER

The brake lever allows switching between “electric drive” mode and “manual pushing” mode.

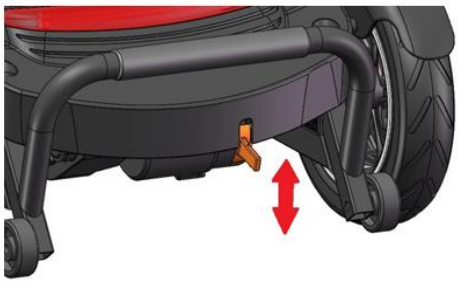
Electric drive: set the lever to position D (locked position) – the vehicle moves using the motor.
Manual mode: set the lever to position N (unlocked position) – allows the vehicle to be moved manually.
If the vehicle does not move and only an audible signal is heard, check whether the brake lever is in position D.
WARNING! Do not leave the brake lever in the unlocked (N) position while driving downhill – the brake will not function.
Note: During use, use only the power supply from the left battery.

BATTERY INSTALLATION AND REMOVAL

Removal: Grip the battery handle. Due to the hook-and-loop fasteners on the bottom of the battery housing, gently move the battery forward and backward and side to side before lifting it up.
Installation: Position the battery housing correctly relative to the socket at the bottom and slide it into place.
WARNING! During transport, always make sure that the AT52339 mobility scooter and its battery are properly secured. The batte ry must be kept in an upright position and must not be transported together with flammable materials.
Note: To prevent power disconnection while driving (e.g. on uneven surfaces), the battery is mounted in a special way. After installation, it must be pressed down firmly so that the hook-and-loop fasteners on the bottom adhere properly.
Note: Due to the hook-and-loop fasteners located on the bottom of the housing, greater force is required for removal and the battery must be pulled upward.
Note: Make sure that the battery cables are connected correctly.

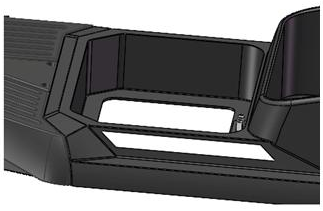FOLDING, UNFOLDING, AND TRANSPORT OF THE SCOOTER

The scooter with removable battery consists of four main parts: frame, battery, seat, and basket. During disassembly, follow the instructions and exercise caution.

Folding procedure:

Folding the scooter

With one hand, gently move the steering column while turning the knob at the same time – this will make it easier to tighten it.

With one hand, lift the rear part of the seat, grasp its center, and lift it upward.


While lifting the seat, the folding lever will automatically be locked by the spring on the pin – you will hear a "click."

Press the seat to make sure that it is properly locked.

Fold the armrest forward until resistance is felt.

WARNING! The folding lever must be properly locked by the pin. Otherwise, there is a risk of falling while seated.

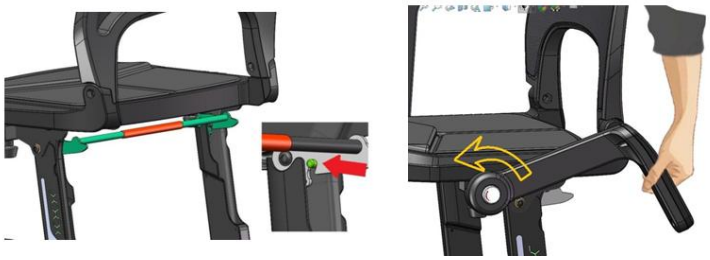
Unfolding procedure:

Perform the steps in reverse order.

With one hand, lift the folding lever, and with the other hand press the seat downward.

With one hand, press the button, and with the other hand rotate the armrest backward.

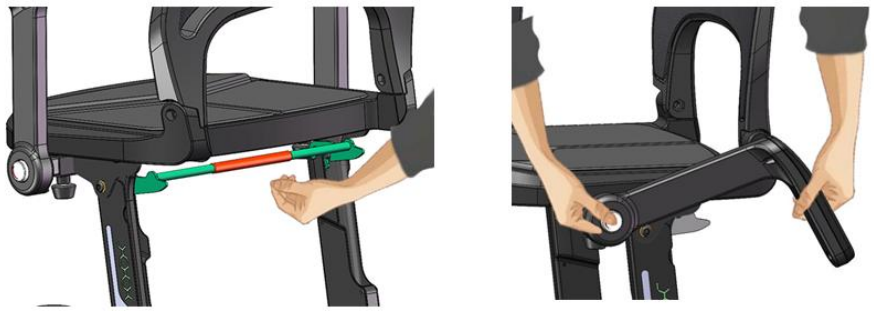Transport and carrying: When carrying the scooter, hold it by the front and rear protective bars (bumpers), as shown below.

CHARGING THE BATTERIES

Insert the charger plug (8) into the charging socket (10) located on the housing.


Connect the charger power cord (9) to the mains socket.

After charging is completed, disconnect the power cord (9) from the socket, and then remove the charger plug (8) from socket (10).

Under ideal storage conditions, fully charged and unused batteries should be recharged every 20 days.

If you plan a longer break from using the AT52339 mobility scooter, it is recommended to charge it for two days and then disconnect the batteries.

If the scooter has not been used for a long time, the batteries should be charged for at least 24 hours before driving.

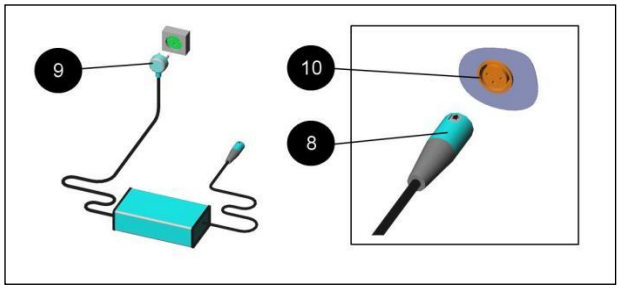Recommended charging cycle

Use the scooter during the day as needed or until the battery level is low (according to the LCD indicator).

At the end of the day, charge the batteries overnight.

There is no need to disconnect the batteries after charging is complete – the charger is automatic and stops charging by itself.

There is no need to charge after every ride unless you do not plan to use the scooter again that day.

As long as the charger plug is in the socket, the electronic system of the AT52339 mobility scooter automatically cuts off the power supply – driving is impossible.

The supplied charger is intended for lithium batteries. Use only the recommended type of charger. If a different one is used, consult the seller.

The battery is equipped with a protection circuit. After it is activated, the battery must be recharged to restore its operation.

NOTE: In the event of a "serious incident" related to the device that directly or indirectly led, could have led, or may lead to any of the following:

a) death of a patient, user, or other person, or

b) temporary or permanent deterioration of the health of a patient, user, or other person, or

c) a serious public health threat

such a "serious incident" must be reported to the manufacturer and to the competent authority of the Member State in which the user or patient resides. In the United Kingdom, the competent authority is the Medicines and Healthcare products Regulatory Agency (MHRA).

NOTE: In the event of pain, allergic reactions, or other disturbing or unclear symptoms related to the use of the medical device, consult a healthcare professional. NOTE: It is prohibited to use the product in any manner other than its intended purpose.

CLEANING, MAINTENANCE AND STORAGE

Cleaning

The AT52339 mobility scooter should be cleaned as needed using a soft, damp cloth and mild detergents. After cleaning, wipe the surface with clean water and then dry thoroughly.

Do not use pressure washers or large amounts of water, as this may damage electrical components. Cleaning should be carried out at temperatures above 0°C.

Maintenance

The scooter requires regular inspections despite low maintenance requirements:

Daily: check the battery charge level, horn operation, and cleanliness of the battery compartment.

Weekly: check the operation of the drive lever, brake, and the stability of wheels and connections.

Monthly: check tire wear and proper driving and braking performance.

Periodically: inspect the condition of cables, connections, and all structural components.

All repairs, adjustments, and battery replacement should be carried out only by qualified service personnel.

Storage

In case of prolonged non-use (e.g., over 4 months), the scooter should be stored in a dry, clean, and secure place, preferably with access to power supply. It is recommended to cover the scooter during storage. The battery should be recharged at least every 20 days. Do not store the scooter at temperatures below 0°C. Low temperatures (below 10°C) may reduce battery performance and capacity. WARNING! Do not allow the battery to freeze – this may cause permanent damage.

DISPOSAL: At the end of its service life, the scooter must be disposed of in accordance with applicable local, national, and EU regulations regarding waste electrical and electronic equipment (WEEE Directive). The device must not be disposed of with unsorted municipal waste. The scooter contains recyclable components such as plastics, electronic components, and batteries. Particular attention should be paid to the battery, which contains substances requiring appropriate handling (including lithium compounds) and should be delivered to specialized collection points. For information on proper disposal methods, contact local authorities or an authorized seller. NOTE: Improper disposal may pose a risk to the environment and human health.

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